



INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)


(51) International Patent Classification ⁷ : C12Q 1/68		A1	(11) International Publication Number: WO 00/36142
			(43) International Publication Date: 22 June 2000 (22.06.00)
(21) International Application Number: PCT/CA99/01177 (22) International Filing Date: 10 December 1999 (10.12.99) (30) Priority Data: 60/111,794 11 December 1998 (11.12.98) US (71) Applicant (for all designated States except US): VISIBLE GENETICS INC. [CA/CA]; 700 Bay Street, P.O. Box 333, Toronto, Ontario M5G 1Z6 (CA). (72) Inventor; and (75) Inventor/Applicant (for US only): SHIPMAN, Robert [CA/CA]; #2-6240 Montevideo Road, Mississauga, Ontario L5N 2N7 (CA). (74) Agent: DEETH WILLIAMS WALL; National Bank Building, Suite 400, 150 York Street, Toronto, Ontario M5H 3S5 (CA).		(81) Designated States: AU, CA, JP, US, European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE). Published <i>With international search report.</i> <i>Before the expiration of the time limit for amending the claims and to be republished in the event of the receipt of amendments.</i>	

(54) Title: METHOD AND KIT FOR THE CHARACTERIZATION OF ANTIBIOTIC-RESISTANCE MUTATIONS IN *MYCOBACTERIUM TUBERCULOSIS*

(57) Abstract

Amplification and cycle sequencing primer sets have been developed for the detection and analysis of antibiotic resistance-associated mutations in defined regions of the *rpoB* (rifampin), *katG* (isoniazid), *oxyR-ahpC* PR (isoniazid), *mabA* (isoniazid), *rpsL/s12* (streptomycin), *16S/rrs* (streptomycin), *embB* (ethambutol), *pncA* (pyrazinamide), *gyrA* (ciprofloxacin) and *23S* (azithromycin) genes of *Mycobacterium tuberculosis*. These primers can be used in a method for detection and characterization of *Mycobacterium tuberculosis* present in a sample. The method includes the steps of obtaining a sputum sample suspected of containing *M. tuberculosis*, performing a first sequencing procedure, with or without prior amplification, on the sample to detect the presence of *M. tuberculosis*, and if present to evaluate the

rpoB, *katG*, *rpsL/s12* and *23S* genes for the presence of antibiotic-resistance inducing mutations; and (c) if *M. tuberculosis* is detected in step (b), performing a second sequencing procedure, with or without prior amplification, on the sample to evaluate the additional genes for the presence of antibiotic-resistance inducing mutations.

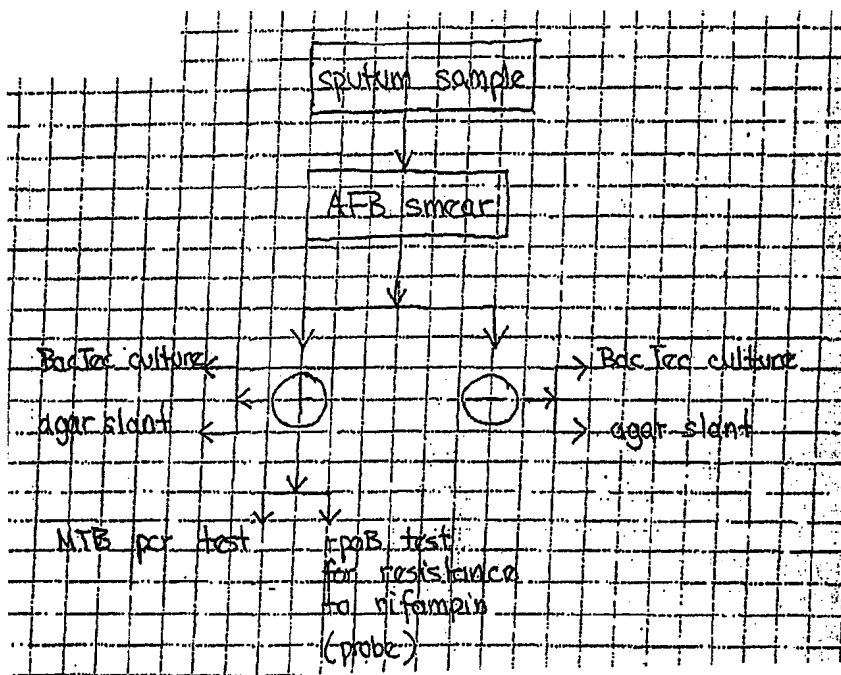

FOR THE PURPOSES OF INFORMATION ONLY

Codes used to identify States party to the PCT on the front pages of pamphlets publishing international applications under the PCT.

AL	Albania	ES	Spain	LS	Lesotho	SI	Slovenia
AM	Armenia	FI	Finland	LT	Lithuania	SK	Slovakia
AT	Austria	FR	France	LU	Luxembourg	SN	Senegal
AU	Australia	GA	Gabon	LV	Latvia	SZ	Swaziland
AZ	Azerbaijan	GB	United Kingdom	MC	Monaco	TD	Chad
BA	Bosnia and Herzegovina	GE	Georgia	MD	Republic of Moldova	TG	Togo
BB	Barbados	GH	Ghana	MG	Madagascar	TJ	Tajikistan
BE	Belgium	GN	Guinea	MK	The former Yugoslav Republic of Macedonia	TM	Turkmenistan
BF	Burkina Faso	GR	Greece	ML	Mali	TR	Turkey
BG	Bulgaria	HU	Hungary	MN	Mongolia	TT	Trinidad and Tobago
BJ	Benin	IE	Ireland	MR	Mauritania	UA	Ukraine
BR	Brazil	IL	Israel	MW	Malawi	UG	Uganda
BY	Belarus	IS	Iceland	MX	Mexico	US	United States of America
CA	Canada	IT	Italy	NE	Niger	UZ	Uzbekistan
CF	Central African Republic	JP	Japan	NL	Netherlands	VN	Viet Nam
CG	Congo	KE	Kenya	NO	Norway	YU	Yugoslavia
CH	Switzerland	KG	Kyrgyzstan	NZ	New Zealand	ZW	Zimbabwe
CI	Côte d'Ivoire	KP	Democratic People's Republic of Korea	PL	Poland		
CM	Cameroon	KR	Republic of Korea	PT	Portugal		
CN	China	KZ	Kazakstan	RO	Romania		
CU	Cuba	LC	Saint Lucia	RU	Russian Federation		
CZ	Czech Republic	LI	Liechtenstein	SD	Sudan		
DE	Germany	LK	Sri Lanka	SE	Sweden		
DK	Denmark	LR	Liberia	SG	Singapore		
EE	Estonia						

- 1 -

**METHOD AND KIT FOR THE CHARACTERIZATION OF
ANTIBIOTIC-RESISTANCE MUTATIONS IN
*MYCOBACTERIUM TUBERCULOSIS***

DESCRIPTION

Field of the Invention

This application relates to a method and kit for the characterization of antibiotic resistance mutations in *Mycobacterium tuberculosis*, and particularly to the evaluation of such mutations in clinical samples.

Background of the Invention

M. tuberculosis can be resistant to all antibiotics that are currently used to treat tuberculosis patients. Antibiotic resistance is due to acquired point mutations in target genes in the genome of *M. tuberculosis*. These point mutations render the organism insensitive to the action of the antibiotic by preventing its uptake or activation, or by altering the antibiotic target. The observed antibiotic resistance in *M. tuberculosis* is not due to an episome-encoded resistance gene transferred from one strain to another and, like other bacteria, is single-step (one point mutation), high-level resistance.

Rapid and accurate detection of antibiotic resistance in *Mycobacterium tuberculosis* in sputum samples would greatly improve both patient treatment and outcome. Presently, analysis of *M. tuberculosis* is carried out on DNA recovered from sputum samples handled according to Standard Infectious Disease/Public Health Laboratory practices. The sputum sample is decontaminated and a cell sediment isolated. This cell sediment is the sample source for all routine procedures used in the detection and isolation of *M. tuberculosis*. Portions of this sample are used in BacTec cultures for selective growth of *M. tuberculosis*, agar plate/agar slant cultures for *M. tuberculosis*, acid-fast bacilli (AFB) smears for mycobacteria detection and molecular biological methods for the detection of *M. tuberculosis* and atypical mycobacteria. (See Fig. 1)

Mycobacterial DNA is prepared directly from the decontaminated sputum cell sediments according to standard procedures and this mycobacterial DNA is used in the various molecular biological detection procedures. The methods presently in use for the detection of *M. tuberculosis* are either PCR-based or probe-based. These tests are used

- 2 -

primarily on AFB smear-positive samples. Since the presence of *M. tuberculosis* has already been established by the AFB smear, these tests are used primarily in a confirmatory capacity as opposed to a diagnostic capacity. Furthermore, these tests provide no information on the potential antibiotic resistance of these *M. tuberculosis* samples.

Below is a list of antibiotics used to treat *M. Tuberculosis* infections. The gene target of the specific antibiotic and regions associated with antibiotic resistance are listed, if known. The references on which the codon assignments are based are listed at the end of the application.

1.	Rifampin	rpoB gene	codon 507-533 ^a
2.	Isoniazid	katG gene	codon 275/315/328 ^b
3.	Isoniazid	mabA gene	unknown ^c
4.	Isoniazid	oxyR-ahpC intergenic region (PR)	nucleotides -48 to +33
5.	Azithromycin	23S rRNA sequence	nucleotide 2568A ^e
6.	Pyrazinamide	pncA gene	codon 47/85 ^f
7.	Ethambutol	embB gene	codon 306 ^g
8.	Streptomycin	rpsL/s12 gene	codon 43/88 ^h
9.	Streptomycin	16S/rrs sequence	nucleotides 491, 512, 516, 513, 903, 904
10.	Ciprofloxacin	gyrA gene	codon 88-95 ^j

Probe-based tests do exist for the determination of rifampin resistance in *M. tuberculosis* (line probe assay-InnoTek), but these probes rely on prior knowledge of antibiotic resistance-associated mutations in the rpoB gene. Mutations outside the region covered by the probe or new mutations not included in the probe cocktail could still confer resistance, but would not be detected using this product in its present form.

Thus, there remains a need for a method for detecting antibiotic-resistance mutations in clinical *M. tuberculosis* sputum samples which is capable of detecting mutations in all of the gene targets which confer antibiotic resistance. It is an object of the present invention to provide such a method. It is a further object of this invention to provide amplification and cycle sequencing primer sets, and kits containing such primer sets, for use in the characterization of antibiotic resistance mutations in *M. tuberculosis*.

Summary of Invention

Amplification and cycle sequencing primer sets have been developed for the detection and analysis of antibiotic resistance-associated mutations in defined regions of the *rpoB* (rifampin), *katG* (isoniazid), *oxyR-ahpC PR* (isoniazid), *mabA* (isoniazid), *rpsL/s12* (streptomycin), *16S/rrs* (streptomycin), *embB* (ethambutol), *pncA* (pyrazinamide), *gyrA* (ciprofloxacin) and *23S* (azithromycin) genes. Using these primer sets and the OPENGENE™ automated DNA sequencing system, a protocol has been developed which permits both the rapid detection of *M. tuberculosis* and the identification of antibiotic resistance-associated mutations in a series of gene targets. The present invention uses a series of tests designed to detect antibiotic resistance-associated mutation in all gene targets for all antibiotics presently in use for the treatment of tuberculosis. The tests are employed in a hierarchical manner on both AFB smear-positive or smear-negative samples to determine both the presence and antibiotic-resistance of *M. tuberculosis* in a given sample. This method permits the simultaneous determination of *M. tuberculosis* presence in a sample and the antibiotic-resistance profile to an entire panel of antibiotics. Standard methods require from 2-6 weeks to culture *M. tuberculosis* and additional time to establish antibiotic resistance. Although DNA sequence-based (genotypic) tests are not intended to replace the traditional culture-based (phenotypic) methods, these tests do represent a rapid, sensitive and accurate protocol which provides clinicians with valuable information regarding antibiotic treatment options within days as opposed to weeks.

Brief Description to the Figures

Fig. 1 shows known testing protocols for *M. tuberculosis*; and

Fig. 2 shows a hierarchical assay scheme for evaluating *M. tuberculosis* type in accordance with the invention.

Detailed Description of the Invention

In accordance with the invention, regions of the genome of *M. tuberculosis* associated with antibiotic resistance are amplified and sequenced using specifically designed amplification and sequencing primers. Various techniques for amplification are known, including the basic PCR amplification techniques described in US Patent No. 4,683,202, which is incorporated herein by reference. Similarly, various techniques for sequencing are known, some of which require prior amplification and some of which do not. Included among known sequencing techniques are those disclosed in US Patents Nos. 5,834,189 and 5,789,168, which are incorporated herein by reference. The primers of the invention can be used in any of these sequencing formats, although the invention is exemplified below using separate amplification and cycle-sequencing steps.

In theory, the selection of primers to amplify and sequence a known region of interest should be straightforward. In fact, however, because of the possibility of primer binding to other sites, complications arising from secondary structure, and other factors which are not fully understood, some primers perform better than others for amplification and sequencing of the same region of interest. The present invention provides primers which have been optimized for the amplification and sequencing of regions associated with each of the ten known types of antibiotic resistance. These primer sets are shown below, along with the sequence of the genes that they are used to analyze. In the gene sequences, the locations of the primers are underlined.

Primers

rpoB (rifampin resistance)

rpoB-F amplification primer, 20-mer, bp2201-2220

5' TAC GGT CGG CGA GCT GAT CC 3'

SEQ. ID NO. 1

rpoB-R amplification primer, 20-mer, bp2611-2592

5' TAC GGC GTT TCG ATG AAC CC 3'

SEQ ID NO. 2

- 5 -

rpoB-5s sequencing primer, 20-mer, bp2201-2220

5' TAC GGT CGG CGA GCT GAT CC 3'

SEQ ID NO. 3

rpoB-3s sequencing primer, 20-mer, bp2611-2592

5' TAC GGC GTT TCG ATG AAC CC 3'

SEQ ID NO. 4

SEQ. ID. NO. 5

2161 aaaccgacga catcgaccac ttccggcaacc gccgcctgcg ~~tacggtcggc gagctgatcc~~
 2221 aaaaccagat ccgggtcggc atgtcgcgga tggagcgggt ggtccgggag cggatgacca
 2281 cccaggacgt ggaggcgatc acaccgcaga cgttgatcaa catccggccg gtggtcgccg
 2341 cgatcaagga gttcttcggc accagccagc tgagccaatt catggaccag aacaaccgcg
 2401 tgtcgggggtt gaccacaaag cgcgcactgt cggcgctggg gcccgggcgt ctgtcacgtg
 2461 agcgtgccgg gctggaggtc cgcgacgtgc acccgtcgca ctacggccgg atgtgcccga
 2521 tcgaaacccc tgagggggccc aacatcggtc tgatcggtc gctgtcggtg tacgcgcggg
 2581 tcaacccggt cggttcacac gaaacgccgt accgcaaggt ggtcgacggc gtgggttagcg

katG (isoniazid resistance)

katG-F amplification primer, 20-mer, bp722-741

5' ATG GGG CTG ATC TAC GTG AA 3'

SEQ ID NO. 6

katG-R amplification primer, 20-mer, bp1250-1231

5' GGT GTT CCA GCC AGC GAC GC 3'

SEQ ID NO. 7

katG-5s sequencing primer, 20-mer, bp722-741

5' ATG GGG CTG ATC TAC GTG AA 3'

SEQ ID NO. 8

katG-3s sequencing primer, 20-mer, bp1250-1231

5' GGT GTT CCA GCC AGC GAC GC 3'

SEQ ID NO. 9

SEQ ID NO. 10

661 gctcggcgat gagcggttaca gcggtaagcg ggatctggag aaccgcgtgg ccgcgggtgca
 721 gatggggctg atctacgtga acccggaggg gccgaacggc aaccgggacc ccatggccgc
 781 ggcggtcgac attcgcgaga cgtttcggcg catggccatg aacgacgtcg aaacagcggc
 841 gctgatcgtc ggcggtcaca ctttcggtaa **gacc**catggc gccggccccg ccgatctggt
 901 cggccccgaa cccgaggctg ctccgctgga gcagatgggc ttgggctgga agagctcgta
 961 tggcaccgga accggtaagg acgcgatcac **cagg**cggcatc gaggtcgat ggacgaacac
 1021 cccgacgaaa **tggg**gacaaca gtttctcga gatcctgtac ggctacgagt gggagctgac

- 6 -

1081 gaagagccct gctggcgctt ggcaatacac cgccaaggac ggcgccggtg ccggcaccat
 1141 cccggaccgcg ttccggcgggc cagggcgctc cccgacgatg ctggccactg acctctcgct
 1201 gcgggtggat ccgatctatg agcggatcac gcgtcgctgg ctggaacacc ccgaggaatt
 1261 ggccgacgag ttccccaagg cctggtacaa gctgatccac cgagacatgg gtcccgttgc

oxyR-aphC intergenic region (PR)

PR-F amplification primer, 20-mer, bp451-470

5' ACC ACT GCT TTG CCG CCA CC 3'

SEQ ID NO. 11

PR-R amplification primer, 20-mer, bp687-668

5' CCG ATG AGA GCG GTG AGC TG 3'

SEQ ID NO. 12

PR-5s sequencing primer, 20-mer, bp451-470

5' ACC ACT GCT TTG CCG CCA CC 3'

SEQ ID NO. 13

PR-3s sequencing primer, 20-mer, bp687-668

5' CCG ATG AGA GCG GTG AGC TG 3'

SEQ ID NO. 14

SEQ ID NO. 15

361 atgccctggg ggtgcaccga gaccggcttc cgaccaccgc tcgccgcaac gtcgactggc
 421 tcatatcgag aatgcttgcg gcactgctga accactgctt tgccgccacc gcggcgaacg
 481 cgcgaagccc ggccacggcc ggctagcacc tcttggcggc gatgccgata aatatggtgt
 541 gatatatcac ctttgctga cagcgacttc acggcacgat ggaatgtcgc aaccaaatgc
 601 attgtccgct ttgatgatga ggagagtcac gccactgcta accattggcg atcaattccc
 661 cgcctaccag ctacccgctc tcatcgccgg tgacctgtcc aaggtcgacg ccaagcagcc
 721 cggcgactac ttcaccacta tcaccagtga cgaacaccca ggcaagtggc ggggtggtgtt

mabA (isoniazid resistance)

mabA-F amplification primer, 20-mer, bp56-75

5' CCT CGC TGC CCA GAA AGG GA 3'

SEQ ID NO. 16

mabA-R amplification primer, 20-mer, bp303-284

5' ATC CCC CGG TTT CCT CCG GT 3'

SEQ ID NO. 17

mabA-5s sequencing primer, 20-mer, bp56-75

5' CCT CGC TGC CCA GAA AGG GA 3'

SEQ ID NO. 18

- 7 -

mabA-3s sequencing primer, 20-mer, bp303-284

5' ATC CCC CGG TTT CCT CCG GT 3'

SEQ ID NO. 19

SEQ ID NO. 20

```
1   agcgcgacat acctgctgcg caattcgtag ggcgtaata caccgcagc cagggcctcg
61  ctgcccagaa agggatccgt catggtcgaa gtgtgctgag tcacaccgac aaacgtcacg
121 agcgtaaccc cagtgcgaaa gttcccgcgc gaaatcgag ccacgttacg ctggtggaca
181 taccgatttc ggcccggccg cggcgagacg ataggttgtc ggggtgactg ccacagccac
241 tgaagggggc aaacccccat tcgtatcccg ttcagtcctg gttaccggag gaaaccgggg
301 gatcgggctg gcgatcgac agcggctggc tgccgacggc cacaaggtgg ccgtcaccca
```

rpsL/s12 (streptomycin resistance)

s12-F amplification primer, 20-mer, bp1-20

5' CGG TAG ATG CCA ACC ATC CA 3'

SEQ ID NO. 21

s12-R amplification primer, 20-mer, bp384-365

5' GCA TCA GCC CTT CTC CTT CT 3'

SEQ ID NO. 22

s12-5s sequencing primer, 20-mer, bp1-20

5' CGG TAG ATG CCA ACC ATC CA 3'

SEQ ID NO. 23

s12-3s sequencing primer, 20-mer, bp384-365

5' GCA TCA GCC CTT CTC CTT CT 3'

SEQ ID NO. 24

SEQ ID NO. 25

```
1   cggtagatgc caaccatcca gcagctggtc cgcaagggtc gtcgggacaa gatcagtaag
61  gtcaagaccg cggctctgaa gggcagcccg cagcgctcgtg gtgtatgcac ccgctgttac
121 accaccactc cgaagaagcc gaactcggcg cttcggaagg ttgcccgcgt gaagttgacg
181 agtcaggctg aggtcacggc gtacattccc ggcgagggcc acaacctgca ggagcactcg
241 atggtgctgg tgcgcggcgg ccgggtgaag gacctgcctg gtgtgcgcta caagatcatc
301 cgcggttcgc tggatacgca ggggtgtcaag aaccgcaaac aggcacgcag ccgttacggc
361 gctaagaagg agaagggctg atgccacgca aggggcccgc gcccaagcgt ccgttggtca
```

16S/rrs (streptomycin resistance)

16S-F amplification primer, 21-mer, bp5-25

5' GGT GAT CTG CCC TGC ACT TCG 3'

SEQ ID NO. 26

16S-R amplification primer, 21-mer, bp147-127

- 8 -

5' CGT CAC CCC ACC AAC AAG CTG 3' SEQ ID NO. 27
16S-5s sequencing primer, 21-mer, bp5-25
5' GGT GAT CTG CCC TGC ACT TCG 3' SEQ ID NO. 28
16S-3s sequencing primer, 21-mer, bp147-127
5' CGT CAC CCC ACC AAC AAG CTG 3' SEQ ID NO. 29

SEQ ID NO. 30

1 cgtgggtgat ctgccctgca cttcgggata agcctgggaa actgggtcta ataccggata
61 ggaccacggg atgcatgtct tgtggtggaa agcgcttttag cggtgtggga tgagcccgcg
121 gcctatcagc ttgttggtgg ggtgacg

embB (ethambutol resistance)

embB-F amplification primer, 21-mer, bp7761-7781
5' CGG CAA GCT GGC GCA CCT TCA 3' SEQ ID NO. 31
embB-R amplification primer, 21-mer, bp8040-8020
5' AGC CAG CAC ACT AGC CCG GCG 3' SEQ ID NO. 32
embB-5s sequencing primer, 21-mer, bp7761-7781
5' CGG CAA GCT GGC GCA CCT TCA 3' SEQ ID NO. 33
embB-3s sequencing primer, 21-mer, bp8040-8020
5' AGC CAG CAC ACT AGC CCG GCG 3' SEQ ID NO. 34

SEQ ID NO. 35

7741 cggcatgcgc cggtgatcc ~~cggcaagctg ggcacattc~~ accctgaccg acgcctgggt
7801 gatattcggc ttctgtctct ggcattgtcat cggcgcgaaat tcgtcggacg acggctacat
7861 cctgggcatg gcccgagtcg ccgaccacgc cggctacatg tccaactatt tccgtctggtt
7921 cggcagcccc gaggatccct tcggctggta ttacaacctg ctggcgctga tgacccatgt
7981 cagcgacgcc agtctgtgga tgcgctgcc agacctggcc ~~gccgggctag tgtgctggct~~

pncA (pyrazinamide resistance)

pncA-F amplification primer, 20-mer, bp1-20
5' ATG CGG GCG TTG ATC ATC GT 3' SEQ ID NO. 36

- 9 -

pncA-F amplification primer, 20-mer, bp561-542

5' TCA GGA GCT GCA AAC CAA CT 3'

SEQ ID NO. 37

pncA-5s sequencing primer, 20-mer, bp1-20

5' ATG CGG GCG TTG ATC ATC GT 3'

SEQ ID NO. 38

pncA-3s sequencing primer, 20-mer, bp561-542

5' TCA GGA GCT GCA AAC CAA CT 3'

SEQ ID NO. 39

SEQ ID. NO. 40

```

1   atgcgggctg tgatcatcgt cgacgtgcag aacgacttct gcgaggggtgg ctcgctggcg
61  gtaaccgggtg gcgccgcgct ggcccgcgcc atcagcgact acctggccga agcggcggag
121 taccatcacg tcgtggcaac caaggacttc cacatcgacc cgggtgacca cttctccggc
181 acaccggact attcctcgtc gtggccaccg cattgcgtca gcggtactcc cggcgcggac
241 ttccatccca gtctggacac gtcggcaatc gaggcggtgt tctacaaggg tgcctacacc
301 ggagcgtaca gcggcttcga aggagtcgac gagaacggca cgccactgct gaattggctg
361 cggcaacgcg gcgtcgatga ggtcgatgtg gtcggtattg ccaccgatca ttgtgtgcgc
421 cagacggccg aggacgcggt acgcaatggc ttggccacca ggggtgctggg ggacctgaca
481 gcgggtgtgt cggccgatac caccgtcgcc gcgctggagg agatgcgcac cgccagcgtc
541 gagttgggtt gcagctcctg a

```

gyrA (fluoroquinilone/ciprofloxacin resistance)

gyrA-F amplification primer, 20-mer, bp2383-2402

5' CAG CTA CAT CGA CTA TGC GA 3'

SEQ ID NO. 41

gyrA-R amplification primer, 20-mer, bp2702-2683

5' GGG CTT CGG TGT ACC TCA TC 3'

SEQ ID NO. 42

gyrA-5s sequencing primer, 20-mer, bp2383-2402

5' CAG CTA CAT CGA CTA TGC GA 3'

SEQ ID NO. 43

gyrA-3s sequencing primer, 20-mer, bp2702-2683

5' GGG CTT CGG TGT ACC TCA TC 3'

SEQ ID NO. 44

SEQ ID NO. 45

```

2341 cgaccggatc gaaccgggtg acatcgagca ggagatgcag cgcagctaca tcgactatgc
2401 gatgagcgtg atcgtcggcc gcgcgctgcc ggaggtgcgc gacgggctca agcccgtgca
2461 tcgccgggtg ctctatgcaa tgttcgattc cggttccgc ccggaccgca gccacgccaa

```

- 10 -

2521 gtcggccccgg tcggttgccg agaccatggg caactaccac ccgcacggcg acgcgtcgat
 2581 ctacgacagc ctggtgcgca tggcccagcc ctggtcgctg cgctacccgc tggtaggacgg
 2641 ccagggcaac ttcggctcgc caggcaatga cccaccggcg gcgatgaggt acaccgaagc
 2701 ccggctgacc ccgttggcga tggagatgct gagggaaatc gacgaggaga cagtcgattt

23S (macrolide/azithromycin resistance)

23S-F amplification primer, 20-mer, bp2444-2463

5' CGA AAT TCC TTG TCG GGT AA 3'

SEQ ID NO. 46

23S-R amplification primer, 20-mer, bp2683-2664

5' GTA TTT CAA CAA CGA CTC CA 3'

SEQ ID NO. 47

23S-5s sequencing primer, 20-mer, bp2444-2463

5' CGA AAT TCC TTG TCG GGT AA 3'

SEQ ID NO. 48

23S-3s sequencing primer, 20-mer, bp2683-2664

5' GTA TTT CAA CAA CGA CTC CA 3'

SEQ ID NO. 49

SEQ ID NO. 50

2401 gccccagtaa acggcggtgg taactataac catcctaagg tagcgaatt ccttgtcggg
 2461 taagttccga cctgcacgaa tggcgtaacg acttcccaac tgtctcaacc atagactcgg
 2521 cgaaattgca ctacgagtaa agatgctcgt tacgcgcggc aggacgaaaa gaccccgggg
 2581 ccttcactac aacttggtat tgggtgttcgg tacggtttgt gtaggatagg tgggagactt
 2641 tgaagcacag acgccagttt gtgtggagtc gttgttgaaa taccactctg atcgtattgg

To facilitate detection of the sequencing products using real-time fluorescence-based electrophoresis apparatus (for example, a Visible Genetics OPENGENTM sequencer), at least one of the sequencing primers is preferably labeled with a fluorescent label. The label is selected for compatibility with the sequencing apparatus employed, and may be, for example, fluorescein or a cyanine dye such as CY5.0 OR CY5.5.

The primers of the invention are suitably packaged in a kit. This kit will contain individually packaged amplification and sequencing primers sets for each resistance gene to be evaluated by the kit. Thus, the kit of the invention includes at least 4 primers (two amplification and two sequencing primers), and preferably includes the primer sets for a plurality of resistance genes, most preferably the primer sets for all ten resistance genes.

- 11 -

The suitable protocol for the utilization of these primer sets in the evaluation of *M. tuberculosis* in clinical samples utilizes PCR amplification, followed by cycle sequencing. DNA for use in the test is obtained from a sample of sputum (100ul-10ml). The sputum sample is processed according to Standard Infectious Disease/Public Health Laboratory practices (Mycobacteriology Bench Manual, Laboratory Services Branch, December 1997, Ontario Ministry of Health). The sputum sample is homogenized, decontaminated and concentrated. Mycobacterial DNA is prepared directly from a portion of the concentrated cell sediment (100-200ul) using standard DNA extraction methods or commercially available kits.

Amplification of the DNA is performed using the amplification primer sets described above. PCR reagents can be prepared for individual reactions, or may be prepared as a master mix which can be used for multiple tests e.g., 10 PCR reactions. Exemplary combinations of reagents are summarized in the following table.

PCR mix		1 PCR	10 PCRs	final conc. / PCR
genomic DNA	(20ng/ul)	1.0ul		20ng
(~0.5fM)				
10X PCR buffer I		2.5ul	25.0ul	1X
2.5mM dNTP mix	(1:1:1:1)	2.5ul	25.0ul	250uM
DMSO		1.3ul	13.0ul	5%
Taq DNA polymerase (1U)		0.2ul	2.0ul	1 unit
molecular grade water		16.5ul		165.0ul
MTB gene primers	(10uM)	1.0ul	10.0ul	10pmol per primer
total volume per PCR		25.0ul		

If the master mix as shown in the column labeled 10 PCRs is utilized, the mastermix contains all the necessary PCR reagents other than the genomic DNA. In this example, 24.0ul of the mastermix is added to a PCR tube, that already contains 1.0ul of genomic DNA, prior to the addition of the mineral oil overlay and placement in the thermocycler.

The genomic DNA preparation utilized must be of sufficient quality and integrity for robust and reproducible PCR. Suitable DNA preparation can be obtained using

- 12 -

the Gentra Puregene™ DNA isolation kit. The kit components are appropriate for the isolation of genomic DNA from blood, fresh or frozen tissue, archival material and paraffin-embedded tissue.

Each primer pair is used to amplify a single gene region under the following conditions:

1.	Denaturation	94°C	5 minutes	1 cycle
2.	Denaturation	94°C	30 seconds	
	Annealing	60°C	30 seconds	35 cycles
	Extension	72°C	60 seconds	
3.	Extension	72°C	5 minutes	1 cycle
4.	Hold	6°C		

The temperature change during the cycles of the step 2 is desirably set to ramp at a rate of 1°C/sec.

After amplification, 2.0ul from the 25.0ul PCR is analysed for purity on a 0.8% agarose gel. Samples displaying single PCR product bands can be used directly for sequence analysis. The yield and purity of the PCR product determines the amount to be used in the subsequent cycle sequencing reaction. Comparable verification of sequencing purity is performed on each of the other amplification products.

Sequence analysis is carried out on the amplified product. The basic procedures and conditions are the same for each region. Accordingly, the invention will be exemplified using the rpoB gene.

For initial sequence analysis of rpoB, the rpoB-5s primer should be used. For confirmatory sequence analysis the rpoB-3s primer should be used. For each template to be sequenced, aliquot 3.0ul of each of the nucleotide termination mixes into four separate tubes marked <A>, <C>, <G> and <T> and store on ice until the sequencing mastermix is prepared.

Cycle sequencing mastermix

rpoB template	2.0ul
10X VGI Sequenace™ buffer	2.5ul

- 13 -

DMSO	3.5ul
2.5uM dye-sequencing primer	2.0ul
PCR grade water	9.0ul
<u>1:10 diluted Thermosequenase</u>	<u>0.5 ul</u>
total volume	22.0ul

Mix the DMSO and other components in the mastermix well by repeated pipetting (5 times) with a micropipette. Store the mastermix on ice until ready to add to the nucleotide termination mixes.

Add 5.0ul of the mastermix to each of the four marked tubes containing the nucleotide termination mixes.

Add 8.0ul lightweight mineral oil to each of the four marked tubes containing the mastermix and nucleotide termination mixes.

Store on ice until ready to load into the thermocycler.

Parameters for cycle sequencing

1.	Denaturation	94°C	5 minutes	1X
2.	Denaturation	94°C	30 seconds	
	Annealing	60°C	30 seconds	35X
	Extension	72°C	60 seconds	
3.	Extension	72°C	5 minutes	1X
4.	Hold	6°C		

The temperature change during the cycles of the step 2 is desirably set to ramp at a rate of 1°C/sec.

At the end of the cycle sequencing reaction add 6.0ul of the Stop Loading Dye directly to each of the four tubes to stop the sequencing reaction. The sequencing samples are heated at 95°C for 2 minutes and then placed on ice before loading 2.0ul (from a total volume of 14ul) on the CLIPPER™ sequencer. The remainder of the sequencing reaction can be stored at -20°C for subsequent use.

The CLIPPER™ sequencer is set-up as described in the *OPENGENE Automated DNA Sequencing System User Manual*. Run parameters for the CLIPPER™ sequencer are 54°C/ 1300volts/ 0.5sec sampling/35min run/50% laser power. The samples loaded included 2 ul each of the forward and reverse sequencing reaction products for the target gene, differentially labeled, for example with CY5.0 and CY5.5 cyanine dye labels. Once the run is completed, the base-called data is analysed by comparison of the test sequence to the rpoB sequence database in GENELIBRARIAN™. This sequence alignment compares the test sequence to the standard control sequence and allows sequence ambiguities to be assessed. Once edited the test sequence can be screened for antibiotic resistance-associated mutations using GENELIBRARIAN™.

Testing for multiple types of antibiotic-resistance mutations can be carried out using a hierarchical assay, as summarized in Fig. 2. At present molecular biological methods for the detection of *M. tuberculosis* are only performed on AFB smear-positive sputum samples. These methods serve as confirmatory tests for the presence of *M. tuberculosis*. In addition to these molecular biological methods, the culture-based procedures for *M. tuberculosis* detection (BacTec liquid culture, agar plate and slant cultures) are performed in parallel. AFB smear-negative sputum samples are processed with only the culture-based detection procedures (Figure 1).

In the present invention both AFB smear-positive and smear-negative sputum samples can be processed using both culture-based and molecular biological methods. A limitation of the AFB stain methodology is its limit of detection. If a sputum sample has a mycobacterial concentration of less than 5000 bacteria/ul the AFB stain will be negative. In addition to this is the observation that the decontamination procedure used to prepare the sputum sample usually kills 10-20% of the mycobacteria present. This would suggest that two-thirds of the AFB smear-negative samples potentially contain mycobacteria. In practice 10-20% of the AFB smear-negative samples are culture-positive for *M. tuberculosis* (Ontario Public Health Laboratory). This level of mycobacteria is easily detected by molecular biological methods and is therefore incorporated in the present invention.

The hierarchy proposed incorporates tests that specifically detect *M. tuberculosis* (rpoB), detect mutations in genes associated with resistance to the "first-line" antibiotics used to treat *M. tuberculosis* infections (rpoB, katG, rpsL/s12, PR, embB, pncA)

and detect other species of mycobacteria (23S) in the absence of *M. tuberculosis* (Figure 2). Group I analyses are performed before both Group II and Group III. Group I analysis will provide information on the antibiotic resistance status to rifampin (rpoB), isoniazid (katG), streptomycin (rpsL/s12) and azithromycin (23S). In addition the rpoB amplification indicates the presence of *M. tuberculosis* and in the absence of rpoB amplification the 23S sequence allows identification of most of the clinically relevant mycobacterial species. Group II analysis provides information on antibiotic resistance mutations in the "second-line" antibiotics used to treat *M. tuberculosis* infections namely, isoniazid (PR), ethambutol (embB), pyrazinamide (pncA) and ciprofloxacin (gyrA). Group III contains gene targets in which mutations associated with antibiotic resistance are infrequently found. This protocol permits specific gene targets to be examined according to the local treatment procedures since the both antibiotics used to treat *M. tuberculosis* infections, and thus the associated antibiotic resistance mutation patterns, vary geographically. As shown in Figure 2 the culture-based methods are performed in parallel. The molecular biological methods would permit the identification of *M. tuberculosis* from both AFB smear-positive and smear-negative sputum samples and further provide information on the antibiotic resistance profile of these samples well in advance of current culture-based methods. This information would be crucial to the initiation of appropriate and effective antibiotic treatment regimens for *M. tuberculosis* infections.

Examples

A pool of DNA samples from antibiotic-sensitive *M. tuberculosis* isolates was obtained from the LCDC, Health and Welfare Canada. Ottawa, Ontario. Wild-type sequence traces, for all gene targets known to harbor mutations in antibiotic-resistant *M. tuberculosis*, were generated.

A panel of DNA samples from five phenotypic streptomycin-resistant *M. tuberculosis* isolates was obtained from the Public Health Laboratory, Ontario Ministry of Health, Toronto, Ontario. These DNA samples were examined for antibiotic resistance-associated mutations in all 10 antibiotic gene targets listed above. Streptomycin resistance-associated mutations were detected in the rpsL/s12 gene in four isolates. Parallel antibiotic resistance-associated mutations in the rpoB (rifampin), katG (isoniazid), PR (isoniazid),

embB (ethambutol), pncA (pyrazinamide) and gyrA (ciprofloxacin) genes were also identified which underscores the importance of examining all the gene targets for first-line antibiotics used in the treatment of *M. tuberculosis*. A summary of the results is shown in Table 1.

The following references are cited herein and are incorporated herein by reference for all states which allow such incorporation.

- ^a DL Williams et al. (1994). Characterisation of rifampin resistance in pathogenic mycobacteria. Antimicrob Agents Chemother 38: 2380-2386.**
- ^b WH Haas et al. (1997). Molecular analysis of katG gene mutations in strains of Mycobacterium tuberculosis complex from Africa. Antimicrob Agents Chemother 41: 1601-1603.**
- ^c S Sreevatsan et al. (1997). Analysis of the oxyR-ahpC region in isoniazid-resistant and -susceptible Mycobacterium tuberculosis complex organisms recovered from diseased humans and animals in diverse localities. Antimicrob Agents Chemother 41: 600-606.**
- ^d A Telenti et al. (1994). Genotypic assessment of isoniazid and rifampin resistance in Mycobacterium tuberculosis: a blind study at the reference laboratory level. Antimicrob Agents Chemother 35: 719-723.**
- ^e C Katsukawa et al. (1997). Characterisation of the rpsL and rrs genes of streptomycin-resistant clinical isolates of Mycobacterium tuberculosis in Japan. J Appl Microbiol 83: 634-640.**
- ^f C Katsukawa et al. (1997). Characterisation of the rpsL and rrs genes of streptomycin-resistant clinical isolates of Mycobacterium tuberculosis in Japan. J Appl Microbiol 83: 634-640.**
- ^g MA Lety et al. (1997). A single point mutation in the embB gene is responsible for resistance to ethambutol in Mycobacterium smegmatis. Antimicrob Agents Chemother 41: 2629-2633.**
- ^h A Scorpio et al. (1997). Characterisation of pncA mutations in pyrazinamide-resistant Mycobacterium tuberculosis. Antimicrob Agents Chemother 41: 540-543.**

- 17 -

- i **C Xu et al. (1996). Fluoroquinilone resistance associated with specific gyrase mutations in clinical isolates of multidrug-resistant *Mycobacterium tuberculosis*. J Infect Disease 174: 1127-1130.**
- J **KA Nash et al. (1995). Genetic basis of macrolide resistance in *Mycobacterium avium* isolated from patients with disseminated disease. Antimicrob Agents Chemother 39: 2625-2630.**

Table 1

gene (antibiotic)	OPH#1 bp/codon/aa	OPH#2 bp/codon/aa	OPH#3 bp/codon/aa	OPH#4 bp/codon/aa	OPH#11 bp/codon/aa
rpoB (rifampin)	cac526tac, His526Tyr	tcg553ttg, Ser553Leu	cac526gac, His526Asp	tcg553ttg, Ser553Leu	wt
katG.1 (isoniazid)	agc513acc, Ser513Thr	agc513acc, Ser513Thr	agc513acc, Ser513Thr	wt	wt
oxyR-ahpC PR (isoniazid)	g541a	wt	wt	wt	g541a
fabG (isoniazid)	wt	wt	wt	wt	wt
rpsL/s12 (streptomycin)	wt	aag43agg, Lys43Arg	aag43agg, Lys43Arg	aag88agg, Lys88Arg	aag43agg, Lys43Arg
16s/rrs (streptomycin)	wt	wt	wt	wt	wt
embB (ethambutol)	wt	gtc292ttc, val292phe	wt	wt	wt
pncA (pyrazinamide)	tcc65lct, Ser65Ser	wt	att133aat, Ile133Asn	wt	tcc65lct, Ser65Ser
gyrA (ciprofloxacin)	agc95acc, Ser95Thr	agc95acc, Ser95Thr	agc95acc, Ser95Thr	agc95acc, Ser95Thr	agc95acc, Ser95Thr
23s (azithromycin)	wt	wt	wt	wt	wt

CLAIMS

1. A method for detection and characterization of *Mycobacterium tuberculosis* present in a sample, comprising the steps of:
 - (a) obtaining a sputum sample suspected of containing *M. tuberculosis*,
 - (b) performing a first sequencing procedure, with or without prior amplification, on the sample to detect the presence of *M. tuberculosis*, and if present to evaluate the rpoB, katG, rpsL/s12 and 23S genes for the presence of antibiotic-resistance inducing mutations; and
 - (c) if *M. tuberculosis* is detected in step (b), performing a second sequencing procedure, with or without prior amplification, on the sample to evaluate the additional genes for the presence of antibiotic-resistance inducing mutations.
2. The method of claim 1, wherein the second sequencing procedure evaluates PR, embB pncA and gyrA genes for the presence of antibiotic-resistance mutations.
3. The method of claim 3, further comprising the step of performing a third sequencing procedure when *M. tuberculosis* was detected in step (b) to evaluate 16S/rrs and mabA genes for the presence of antibiotic-resistance mutations.
4. The method of any of claims 1 to 3, wherein the first sequencing procedure for rpoB is performed using amplification primers as set forth in Seq. ID Nos. 1 and 2 and sequencing primers as set forth in Seq. ID. Nos. 3 and 4.
5. The method of any of claims 1 to 4, wherein the first sequencing procedure for katG is performed using amplification primers as set forth in Seq. ID Nos. 6 and 7 and sequencing primers as set forth in Seq. ID. Nos. 8 and 9.
6. The method of any of claims 1 to 5, wherein the first sequencing procedure for rpsL/s12 is performed using amplification primers as set forth in Seq. ID Nos. 21 and 22 and sequencing primers as set forth in Seq. ID. Nos. 23 and 24.

- 19 -

7. The method of any of claims 1 to 6, wherein the second sequencing procedure for 23S is performed using amplification primers as set forth in Seq. ID Nos. 46 and 47 and sequencing primers as set forth in Seq. ID. Nos. 48 and 49.

8. The method of any of claims 1 to 7, wherein the second sequencing procedure for PR is performed using amplification primers as set forth in Seq. ID Nos. 11 and 12 and sequencing primers as set forth in Seq. ID. Nos. 13 and 14.

9. The method of any of claims 1 to 8, wherein the second sequencing procedure for pncA is performed using amplification primers as set forth in Seq. ID Nos. 36 and 37 and sequencing primers as set forth in Seq. ID. Nos. 38 and 39.

10. The method of any of claims 1 to 9, wherein the second sequencing procedure for embB is performed using amplification primers as set forth in Seq. ID Nos. 31 and 32 and sequencing primers as set forth in Seq. ID. Nos. 33 and 34.

11. The method of any of claims 1 to 10, wherein the second sequencing procedure for gyrA is performed using amplification primers as set forth in Seq. ID Nos. 41 and 42 and sequencing primers as set forth in Seq. ID. Nos. 43 and 44.

12. The method of any of claims 2 to 11, wherein the third sequencing procedure for 16S/rrs is performed using amplification primers as set forth in Seq. ID Nos. 26 and 27 and sequencing primers as set forth in Seq. ID. Nos. 28 and 29.


13. The method of any of claims 2 to 12, wherein the third sequencing procedure for mabA is performed using amplification primers as set forth in Seq. ID Nos. 16 and 17 and sequencing primers as set forth in Seq. ID. Nos. 18 and 19.

14. A kit for evaluation of antibiotic-resistance mutations in a sample of *Mycobacterium tuberculosis*, comprising one or more pairs of amplification primers and one or more matched pairs of sequencing primers for amplification and sequencing regions within

the genome of *M. tuberculosis*, characterized in that the amplification and sequencing primer pairs are selected from among:


- (a) amplification primers of Seq. ID Nos. 1 and 2 in combination and sequencing primers of Seq. ID Nos. 3 and 4;
- (b) amplification primers of Seq. ID Nos. 6 and 7 in combination and sequencing primers of Seq. ID Nos. 8 and 9;
- (c) amplification primers of Seq. ID Nos. 11 and 12 in combination and sequencing primers of Seq. ID Nos. 13 and 14;
- (d) amplification primers of Seq. ID Nos. 16 and 17 in combination and sequencing primers of Seq. ID Nos. 18 and 19;
- (e) amplification primers of Seq. ID Nos. 21 and 22 in combination and sequencing primers of Seq. ID Nos. 23 and 24;
- (f) amplification primers of Seq. ID Nos. 26 and 27 in combination and sequencing primers of Seq. ID Nos. 28 and 29;
- (g) amplification primers of Seq. ID Nos. 31 and 32 in combination and sequencing primers of Seq. ID Nos. 33 and 34;
- (h) amplification primers of Seq. ID Nos. 36 and 37 in combination and sequencing primers of Seq. ID Nos. 38 and 39;
- (i) amplification primers of Seq. ID Nos. 41 and 42 in combination and sequencing primers of Seq. ID Nos. 43 and 44; and
- (j) amplification primers of Seq. ID Nos. 46 and 47 in combination and sequencing primers of Seq. ID Nos. 48 and 49.

Figure 1

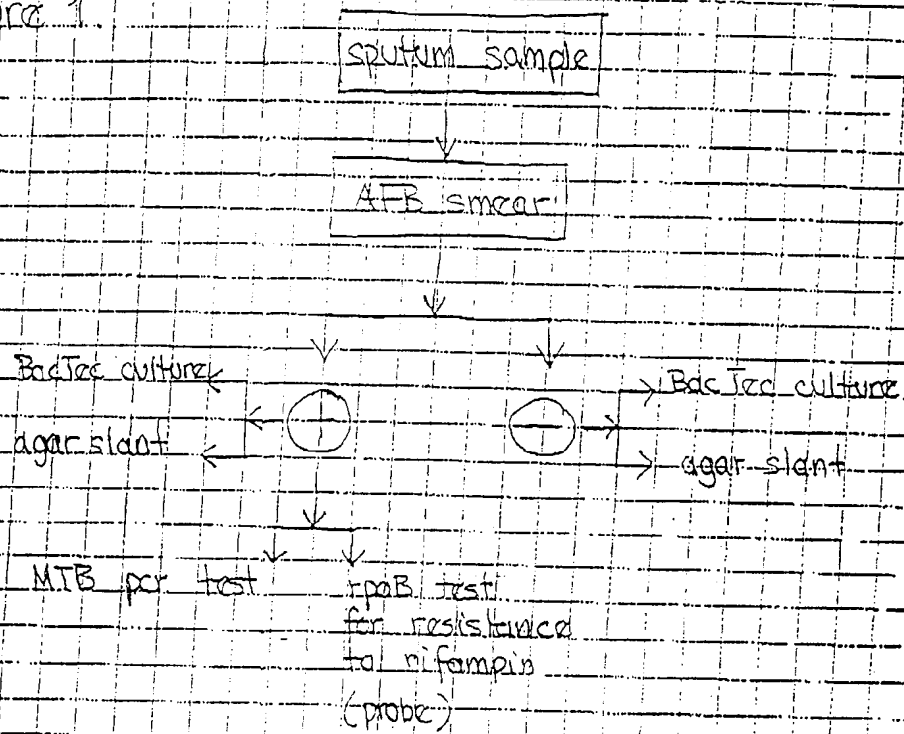

2/2

Figure 2

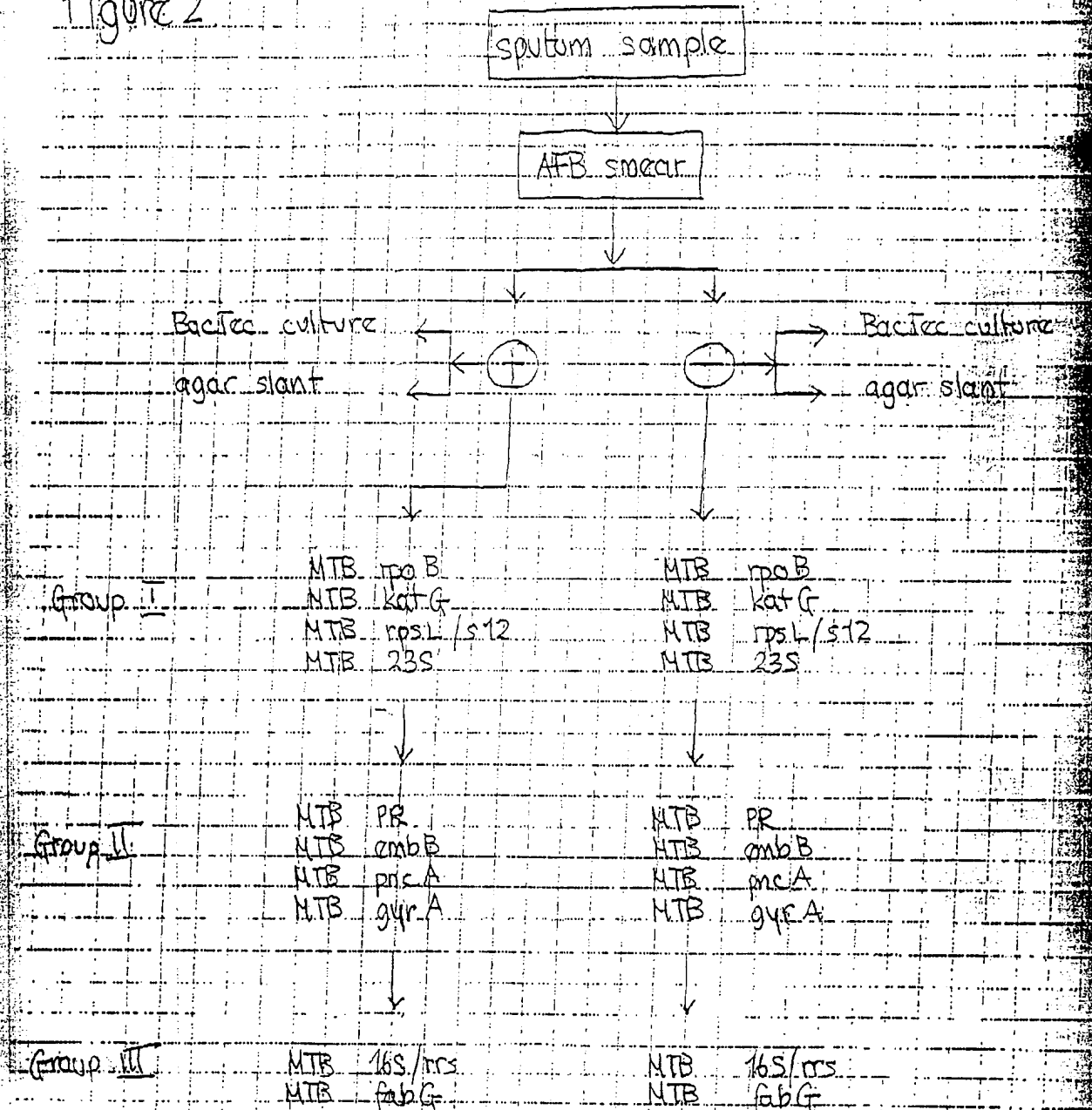

SEQUENCE LISTING¹

<110> Visible Genetics Inc.
Shipman, Robert

<120> Method and Kit for the Characterization of
Antibiotic-Resistance Mutations in Mycobacterium
tuberculosis

<130> VGEN.P-055-WO

<140>

<141>

<150> 60/111,794

<151> 1998-12-11

<160> 50

<170> PatentIn Ver. 2.1

<210> 1

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> rpoB-F amplification primer

<400> 1

tacggtcggc gagctgatcc 20

<210> 2

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> rpoB-R amplification primer

<400> 2

tacggcggtt cgatgaacct 20

<210> 3
<211> 20
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> rpoB-5s sequencing primer

<400> 3
tacggtcggc gagctgatcc 20

<210> 4
<211> 20
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> rpoB-3s sequencing primer

<400> 4
tacggcgttt cgatgaacct 20

<210> 5
<211> 480
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> rpoB (rifampin resistance)

<400> 5
aaaccgacga catcgaccac ttgggcaacc gccgcctgcg tacggtcggc gagctgatcc 60
aaaaccagat ccgggtcggc atgtcgcgga tggagcgggt ggtccgggag cggatgacca 120
cccaggacgt ggaggcgatc acaccgcaga cgttgatcaa catccggccg gtggtcgccg 180
cgatcaagga gttcttcggc accagccagc tgagccaatt catggaccag aacaaccgc 240
tgtcggggtt gaccacaag cgccgactgt cggcgctggg gcccggcggg ctgtcacgtg 300
agcgtgccgg gctggaggtc cgcgacgtgc accgctgca ctacggccgg atgtgcccga 360
tcgaaacccc tgagggggccc aacatcggtc tgatcggtc gctgtcggtg tacgcgcggg 420
tcaaccgtt cgggttcacg gaaacgccgt accgcaaggt ggtcgacggc gtggttagcg 480

<210> 6
<211> 20

<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> katG-F amplification primer

<400> 6
atggggctga tctacgtgaa 20

<210> 7
<211> 20
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> katG-R amplification primer

<400> 7
ggtgttccag ccagcgacgc 20

<210> 8
<211> 20
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> katG-5s sequencing primer

<400> 8
atggggctga tctacgtgaa 20

<210> 9
<211> 20
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> katG-3s sequencing primer

<400> 9
ggtgttccag ccagcgacgc 20

<210> 10
<211> 660
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> katG (isoniazid resistance)

<400> 10
gctcggcgat gagegttaca gcggttaagcg ggatctggag aaccgctgg ccgcggtgca 60
gatggggctg atctacgtga acccggaggg gccgaacggc aaccgggacc ccatggccgc 120
ggcggtcgac attcgcgaga cgttcggcg catggccatg aacgacgtcg aaacagcggc 180
gctgatcgtc ggcggtcaca cttcggtaa gacctatggc gccggcccgg ccgatctggt 240
cggccccgaa cccgaggctg ctccgctgga gcagatgggc ttgggctgga agagctcgta 300
tggcaccgga accggttaagg acgcgatcac cagcggcatc gaggtcgtat ggacgaacac 360
cccgcagaaa tgggacaaca gtttcctcga gatcctgtac ggctacgagt gggagctgac 420
gaagagccct gctggcgctt ggcaatacac cgccaaggac ggcgccggtg ccggcaccat 480
cccggacccg ttcggcgggc cagggcgctc cccgacgatg ctggccactg acctctcgct 540
gcgggtggat ccgatctatg agcggatcac gcgtcgtgg ctggaacacc ccgaggaatt 600
ggccgacgag ttgccaagg cctgtgtaca gctgatccac cgagacatgg gtcccgttc 660

<210> 11
<211> 20
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> PR-F amplification primer

<400> 11
accactgctt tgccgccacc 20

<210> 12
<211> 20
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> PR-R amplification primer

<400> 12
ccgatgagag cggtagctg 20

<210> 13
<211> 20
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> PR-5s sequencing primer

<400> 13
accactgctt tgccgccacc 20

<210> 14
<211> 20
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> PR-3s sequencing primer

<400> 14
ccgatgagag cggtagctg 20

<210> 15
<211> 420
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> oxyR-ahpC intergenic region (PR)

<400> 15
atgccctggg ggtgcaccga gaccggttc cgaccaccgc tcgccgcaac gtcgactggc 60
tcatatcgag aatgcttgcg gcactgctga accactgctt tgccgccacc gcggcgaacg 120
cgcgaagccc ggccacggcc ggctagcacc tcttgccggc gatgccgata aatatggtgt 180
gatatatcac ctttgctga cagcgacttc acggcacgat ggaatgtcgc aaccaaattgc 240
attgtccgct ttgatgatga ggagagtcac gccactgcta accattggcg atcaattccc 300
cgcctaccag ctaccgctc tcacggcgg tgacctgtcc aaggtcgacg ccaagcagcc 360
cggcgactac ttcaccacta tcaccagtga cgaacacca ggcaagtggc gggtggtgtt 420

<210> 16
<211> 20
<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> fabG-F amplification primer

<400> 16

cctcgctgcc cagaaaggga

20

<210> 17

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> fabG-R amplification primer

<400> 17

atcccccggt ttctccggt

20

<210> 18

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> fabG-5s sequencing primer

<400> 18

cctcgctgcc cagaaaggga

20

<210> 19

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> fabG-3s sequencing primer

<400> 19

atcccccggt ttctccggt

20

<210> 20
<211> 360
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> fabG (isoniazid resistance)

<400> 20
agcgcgacat acctgctgcg caattcgtag ggcgtcaata caccgcgagc cagggcctcg 60
ctgcccagaa agggatccgt catggtcgaa gtgtgctgag tcacaccgac aaacgtcacg 120
agcgtaaccc cagtgcgaaa gtcccgccg gaaatcgag ccacgttacg ctcgtggaca 180
taccgatttc ggccggcgcg cggcgagacg ataggtgtc ggggtgactg ccacagccac 240
tgaaggggcc aaacccccat tcgtatccc ttcagtctg gttaccggag gaaaccgggg 300
gatcgggctg gcgatgcac agcggctggc tgccgacggc cacaaggtgg ccgtcaccca 360

<210> 21
<211> 20
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> s12-F amplification primer

<400> 21
cggtagatgc caaccatcca 20

<210> 22
<211> 20
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> s12-R amplification primer

<400> 22
gcatcagccc ttctcttct 20

<210> 23
<211> 20
<212> DNA
<213> Mycobacterium tuberculosis

<220>

<223> s12-5s sequencing primer

<400> 23

cggtagatgc caaccatcca

20

<210> 24

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> s12-3s sequencing primer

<400> 24

gcatcagccc ttctccttct

20

<210> 25

<211> 420

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> rpsL/s12 (streptomycin resistance)

<400> 25

cggtagatgc caaccatcca gcagctggtc cgcaagggtc gtcgggacaa gatcagtaag 60
gtcaagaccg cggctctgaa gggcagcccg cagcgtcgtg gtgtatgcac ccgcgtgtac 120
accaccactc cgaagaagcc gaactcggcg ctctggaagg ttgccgcgt gaagttgacg 180
agtcaggctg aggtcacggc gtacattccc ggcgagggcc acaacctgca ggagcactcg 240
atggtgctgg tgcgcggcgg ccgggtgaag gacctgcctg gtgtgcgcta caagatcatc 300
cgcggttcgc tggatacgca ggggtgcaag aaccgcaaac aggcacgcag ccgttacggc 360
gctaagaagg agaagggtg atgccacgca aggggccccgc gccaagcgt ccgttggtca 420

<210> 26

<211> 21

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> 16S-F amplification primer

<400> 26	
ggtgatctgc cctgcacttc g	21
<210> 27	
<211> 21	
<212> DNA	
<213> Mycobacterium tuberculosis	
<220>	
<223> 16S-R amplification primer	
<400> 27	
cgtcacccca ccaacaagct g	21
<210> 28	
<211> 21	
<212> DNA	
<213> Mycobacterium tuberculosis	
<220>	
<223> 16S-5s sequencing primer	
<400> 28	
ggtgatctgc cctgcacttc g	21
<210> 29	
<211> 21	
<212> DNA	
<213> Mycobacterium tuberculosis	
<220>	
<223> 16S-3s sequencing primer	
<400> 29	
cgtcacccca ccaacaagct g	21
<210> 30	
<211> 147	
<212> DNA	
<213> Mycobacterium tuberculosis	

<220>

<223> 16S/rrs (streptomycin resistance)

<400> 30

cggtgggtgat ctgccctgca ctccgggata agcctgggaa actgggtcta ataccggata 60
ggaccacggg atgcatgtct tgtggtggaa agcgcttag cgggtgtggga tgagcccgcg 120
gcctatcagc ttgttggtgg ggtgacg 147

<210> 31

<211> 21

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> embB-F amplification primer

<400> 31

cggcaagctg gcgcaccttc a 21

<210> 32

<211> 21

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> embB-R amplification primer

<400> 32

agccagcaca ctagcccggc g 21

<210> 33

<211> 21

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> embB-5s sequencing primer

<400> 33

cggcaagctg gcgcaccttc a 21

<210> 34
<211> 21
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> embB-3s sequencing primer

<400> 34
agccagcaca ctagcccggc g 21

<210> 35
<211> 300
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> embB (ethambutol resistance)

<400> 35
cggcatgcgc cggctgattc cggcaagctg ggcaccttc accctgaccg acgccgtggt 60
gatattcggc ttctgtctct ggcatgtcat cggcgcgaat tcgtcggacg acggctacat 120
cctgggcatg gcccagatcg cggaccacgc cggctacatg tccaactatt tccgctgggt 180
cggcagcccg gaggatccct tcggctggta ttacaacctg ctggcgcgtga tgacctatgt 240
cagcgacgcc agtctgtgga tgcgcctgcc agacctggcc gccgggctag tgtgctggct 300

<210> 36
<211> 20
<212> DNA
<213> Mycobacterium tuberculosis

<220>
<223> pncA-F amplification primer

<400> 36
atgcgggcgt tgatcatcgt 20

<210> 37
<211> 20
<212> DNA
<213> Mycobacterium tuberculosis

<220>

<223> pncA-F amplification primer

<400> 37

tcaggagctg caaaccaact

20

<210> 38

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> pncA-5s sequencing primer

<400> 38

atgcgggcgt tgatcatcgt

20

<210> 39

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> pncA-3s sequencing primer

<400> 39

tcaggagctg caaaccaact

20

<210> 40

<211> 561

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> pncA (pyrazinamide resistance)

<400> 40

atgcgggcgt tgatcatcgt cgacgtgcag aacgacttct gcgagggtgg ctgctggcg 60
gtaaccgggtg ggcgcgcgt ggcccgcgc atcagcgact acctggccga agcggcggac 120
taccatcacg tegtggcaac caaggacttc cacatcgacc cgggtgacca cttctccgc 180
acaccggact attcctcgtc gtggccaccg cattgcgtca gcgtactcc cggcgcggac 240
ttccatcca gtctggacac gtcggcaatc gaggcggtgt tctacaaggg tgcctacacc 300

ggagcgtaca gcggcttcga aggagtcgac gagaacggca cgccactgct gaattggctg 360
cggcaacgcg gcgtcgatga ggtcgatgtg gtcggtattg ccaccgatca ttgtgtgcgc 420
cagacggccg aggacgcggt acgcaatggc ttggccacca ggggtctggt ggacctgaca 480
gcgggtgtgt cggccgatac caccgtcgcc gcgctggagg agatgcgcac cgccagcgtc 540
gagttggtt gcagctcctg a 561

<210> 41

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> gyrA-F amplification primer

<400> 41

cagctacatc gactatgcga 20

<210> 42

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> gyrA-R amplification primer

<400> 42

gggcttcggt gtacctcatc 20

<210> 43

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> gyrA-5s sequencing primer

<400> 43

cagctacatc gactatgcga 20

<210> 44

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> gyrA-3s sequencing primer

<400> 44

gggcttcggt gtacctcatc

20

<210> 45

<211> 420

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> gyrA (fluoroquinilone/ciprofloxacin resistance)

<400> 45

cgaccggatc gaaccggtg acatcgagca ggagatgcag cgcagctaca tcgactatgc 60
gatgagcgtg atcgtcggcc gcgcgctgcc ggaggtgcgc gacgggctca agcccggtca 120
tcgccgggtg ctctatgcaa tgttcgattc cggcttcgc ccggaccgca gccacgcaa 180
gtcggcccggtc tcggttgccg agaccatggg caactaccac ccgcacggcg acgcgtcgat 240
ctacgacagc ctggtgcgca tggcccagcc ctggtcgctg cgctacccgc tgggtggacgg 300
ccaggggcaac ttcggctcgc caggcaatga cccaccggcg gcgatgaggt acaccgaagc 360
ccggctgacc ccgttggcga tggagatgct gagggaaatc gacgaggaga cagtcgattt 420

<210> 46

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> 23S-F amplification primer

<400> 46

cgaaattcct tgcgggtaa

20

<210> 47

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> 23S-R amplification primer

<400> 47

gtatttcaac aacgactcca

20

<210> 48

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> 23S-5s sequencing primer

<400> 48

cgaaattcct tgcgggtaa

20

<210> 49

<211> 20

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> 23S-3s sequencing primer

<400> 49

gtatttcaac aacgactcca

20

<210> 50

<211> 300

<212> DNA

<213> Mycobacterium tuberculosis

<220>

<223> 23S (macrolide/azithromycin resistance)

<400> 50

gccccagtaa acggcggtgg taactataac catcctaagg tagcgaaatt ccttgctggg 60
taagttccga cctgcacgaa tggcgtaacg acttccaac tgtctcaacc atagactcgg 120
cgaaattgca ctacagtaaa agatgctcgt tacgcgcggc aggacgaaaa gaccccgga 180
ccttactac aacttggtat tgggttcgg tacggttgt gtaggatagg tgggagactt 240
tgaagcacag acgccagtti gtgtggagtc gttgtgaaa taccactctg atcgattgg 300

INTERNATIONAL SEARCH REPORT

Int. .tional Application No

PCT/CA 99/01177

A. CLASSIFICATION OF SUBJECT MATTER

IPC 7 C12Q1/68

According to International Patent Classification (IPC) or to both national classification and IPC

B. FIELDS SEARCHED

Minimum documentation searched (classification system followed by classification symbols)

IPC 7 C12Q

Documentation searched other than minimum documentation to the extent that such documents are included in the fields searched

Electronic data base consulted during the international search (name of data base and, where practical, search terms used)

C. DOCUMENTS CONSIDERED TO BE RELEVANT

Category °	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
X	KAPUR V ET AL: "Application of automated DNA sequence analysis for mycobacterium species identification and detection of mutations associated with antibiotic resistance in Mycobacterium tuberculosis" ABSTRACTS OF THE INTERSCIENCE CONFERENCE ON ANTIMICROBIAL AGENTS AND CHEMOTHERAPY, vol. 34, 1994, page 163 XP000901974 see abstract D71 ---	1-3
X	SUZUKI Y ET AL: "Detection of kanamycin-resistant Mycobacterium tuberculosis by identifying mutations in the 16SrRNA gene" JOURNAL OF CLINICAL MICROBIOLOGY, vol. 36, no. 5, May 1998 (1998-05), pages 1220-5, XP000901934 the whole document --- -/--	1-3

☒ Further documents are listed in the continuation of box C.

☒ Patent family members are listed in annex.

° Special categories of cited documents :

"A" document defining the general state of the art which is not considered to be of particular relevance

"E" earlier document but published on or after the international filing date

"L" document which may throw doubts on priority claim(s) or which is cited to establish the publication date of another citation or other special reason (as specified)

"O" document referring to an oral disclosure, use, exhibition or other means

"P" document published prior to the international filing date but later than the priority date claimed

"T" later document published after the international filing date or priority date and not in conflict with the application but cited to understand the principle or theory underlying the invention

"X" document of particular relevance; the claimed invention cannot be considered novel or cannot be considered to involve an inventive step when the document is taken alone

"Y" document of particular relevance; the claimed invention cannot be considered to involve an inventive step when the document is combined with one or more other such documents, such combination being obvious to a person skilled in the art.

"&" document member of the same patent family

Date of the actual completion of the international search

13 April 2000

Date of mailing of the international search report

28/04/2000

Name and mailing address of the ISA

European Patent Office, P.B. 5818 Patentlaan 2
NL - 2280 HV Rijswijk
Tel. (+31-70) 340-2040, Tx. 31 651 epo nl,
Fax: (+31-70) 340-3016

Authorized officer

Osborne, H

INTERNATIONAL SEARCH REPORT

International Application No

PCT/CA 99/01177

C.(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT

Category	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
Y	HONORE N ET AL: "Streptomycin resistance in mycobacteria" ANTIMICROBIAL AGENTS AND CHEMOTHERAPY, vol. 38, no. 2, February 1994 (1994-02), pages 238-42, XP000901931 page 239, paragraph 2 ---	1-3
Y	SCORPIO A ET AL: "Characterization of pncA mutations in pyrazinamide-resistant Mycobacterium tuberculosis" ANTIMICROBIAL AGENTS AND CHEMOTHERAPY, vol. 41, no. 3, March 1997 (1997-03), pages 540-543, XP000901990 page 540 -page 542, paragraph 4 ---	1-3
Y	ALANGADEN GJ ET AL: "Mechanism of resistance to amikacin and kanamycin in Mycobacterium tuberculosis" ANTIMICROBIAL AGENTS AND CHEMOTHERAPY, vol. 42, no. 5, May 1998 (1998-05), pages 1295-97, XP000901991 the whole document ---	1-3
Y	HEYM B ET AL: "IMPLICATIONS OF MULTIDRUG RESISTANCE FOR THE FUTURE OF SHORT-COURSE CHEMOTHERAPY OF TUBERCULOSIS: A MOLECULAR STUDY" LANCET THE, GB, LANCET LIMITED. LONDON, vol. 344, no. 8918, 30 July 1994 (1994-07-30), pages 293-298, XP002039609 ISSN: 0140-6736 the whole document ---	1-3
Y	WO 97 23650 A (DUNN JAMES M ; LEUSHNER JAMES (CA); STEVENS JOHN K (CA); VISIBLE GE) 3 July 1997 (1997-07-03) page 10, paragraph 1 -page 11, paragraph 4; example 7 ---	1-3
A	WO 95 33074 A (MAYO FOUNDATION ; HOFFMANN LA ROCHE (US)) 7 December 1995 (1995-12-07) page 3, paragraph 3 -page 6, paragraph 3 ---	1-3
A	WO 95 33851 A (INNOGENETICS NV ; BEENHOUWER HANS DE (BE); PORTAELS FRANCOISE (BE);) 14 December 1995 (1995-12-14) claims 2,22 --- -/--	1-3

INTERNATIONAL SEARCH REPORT

Int. l. Application No

PCT/CA 99/01177

C.(Continuation) DOCUMENTS CONSIDERED TO BE RELEVANT

Category	Citation of document, with indication, where appropriate, of the relevant passages	Relevant to claim No.
P, X	<p>NUESCA D ET AL: "RAPID DETECTION OF ANTIBIOTIC RESISTANCE-ASSOCIATED MUTATIONS IN 10 GENE TARGETS IN MYCOBACTERIUM TUBERCULOSIS USING THE OPENGENE(R) SYSTEM"</p> <p>ABSTRACTS OF THE GENERAL MEETING OF THE AMERICAN SOCIETY FOR MICROBIOLOGY, vol. 99, 30 May 1999 (1999-05-30) - 3 June 1999 (1999-06-03), page 636 XP000891874 see abstract U-13</p> <p>-----</p>	1-3

INTERNATIONAL SEARCH REPORT

Information on patent family members

International Application No

PCT/CA 99/01177

Patent document cited in search report	Publication date	Patent family member(s)	Publication date
WO 9723650 A	03-07-1997	US 5834189 A	10-11-1998
		AU 1426297 A	17-07-1997
		CA 2239896 A	03-07-1997
		EP 0870059 A	14-10-1998
		WO 9741257 A	06-11-1997
		US 5888736 A	30-03-1999
		US 5981186 A	09-11-1999
		US 6007983 A	28-12-1999
WO 9533074 A	07-12-1995	US 5643723 A	01-07-1997
		AU 2605795 A	21-12-1995
WO 9533851 A	14-12-1995	AU 703947 B	01-04-1999
		AU 2789695 A	04-01-1996
		BR 9507960 A	02-09-1997
		CZ 9603612 A	16-07-1997
		EP 0771360 A	07-05-1997
		JP 10500857 T	27-01-1998